


## Restore factory settings AvaSpec


Start the program by double clicking the executable 'restore\_eeprom.exe'

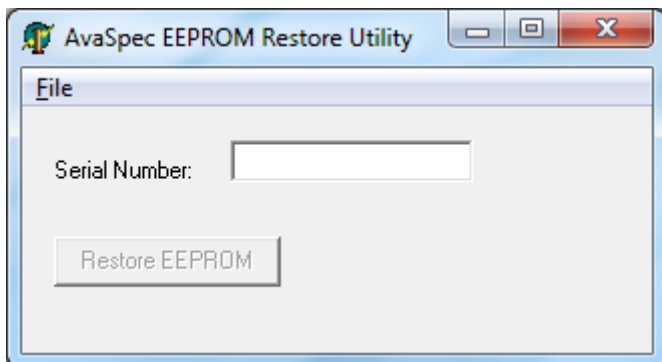

Provide the serial number of the spectrometer that needs to be restored in the input box "Serial Number".

Next load the restore .bak file that is shipped with this tool, with 'File', 'Load .bak file', select the file in the dialog.

Please load only .bak files that were distributed with this utility, otherwise corrupt data might be loaded into the EEPROM.

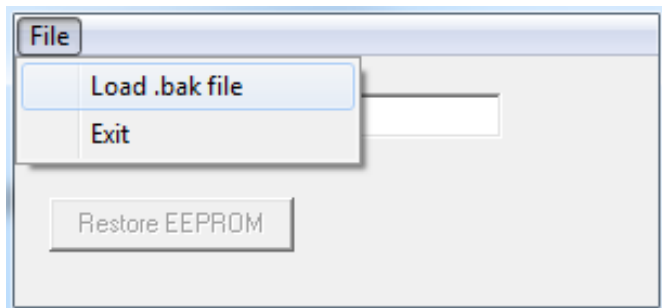

Restore the EEPROM by clicking the button "Restore EEPROM".

The program will report success or failure. On failure please contact Avantes through [support@avantes.com](mailto:support@avantes.com) and mention that you used the restore eeprom program and it failed with a possible error message.

If the program could not find the spectrometer associated with the .bak file, an error message will also be displayed.

Exit the program with 'File', 'Exit', or by clicking the close icon in the upper right corner.

After restoring the EEPROM, please restore the 'Factory Settings' in AvaSoft, by selecting 'Setup', 'Wavelength Calibration Coefficients', clicking the 'Restore Factory Settings' button and clicking the 'OK' button.